



---

**Document Version:** v1.0


**Date of Issue:** June 2025

**Product Model:** Heeve® Aluminium Curved Non-Folding Loading Ramp

**Capacity:** 680kg

**Brand:** Heeve®

**Distributor:** Ramp Champ®

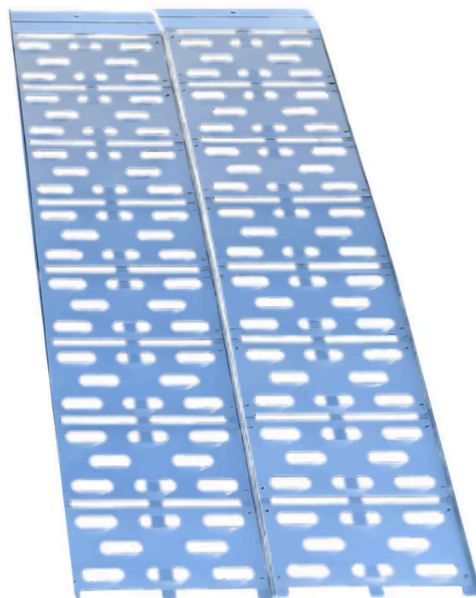

# TABLE OF CONTENTS

- [1. Introduction](#)
- [2. Product Specifications](#)
  - [Gradient Guidelines](#)
- [3. Components & Setup](#)
  - [Setup Instructions:](#)
- [4. Safe Usage Instructions](#)
- [5. Loading & Unloading](#)
  - [5.1 Loading](#)
  - [5.2 Unloading](#)
- [6. Maintenance & Care](#)
- [7. Safety Warnings](#)
- [8. Warranty & Support](#)

# Heeve® Aluminium Curved Non-Folding Loading Ramp

## 1. Introduction

Thank you for purchasing the Heeve® Aluminium Curved Non-Folding Loading Ramp. This access ramp is designed to enable safe mobility for users with wheelchairs, scooters, and/or other equipment. Please read and follow all instructions in this manual carefully before using it.

## 2. Product Specifications

| Heeve® Aluminium Curved Non-Folding Loading Ramp |                                                                                        |
|--------------------------------------------------|----------------------------------------------------------------------------------------|
| Ramp Type                                        | Aluminium Curved Non-Folding Loading Ramp                                              |
| Total Lengths Available                          | 900mm to 1800mm                                                                        |
| Overall Ramp Width                               | 280mm                                                                                  |
| Maximum Capacity                                 | 680kg                                                                                  |
| Ramp Material                                    | Marine-grade aluminium                                                                 |
| Surface                                          | Perforated surface                                                                     |
| Application                                      | Suitable for manual or powered wheelchairs, scooters, trolleys, motorcycles, and more. |

### Gradient Guidelines

- Occupied loading height must be maintained at a maximum of 1:6 gradient or 9.5 degrees.
- Unoccupied loading height should be at a maximum of 1:4 gradient or 14.5 degrees.

| SKU   | Ramp Length | Occupied Max Loading Height | Unoccupied Max Loading Height |
|-------|-------------|-----------------------------|-------------------------------|
| H0022 | 900mm       | 150mm                       | 350mm                         |
| H0023 | 1200mm      | 200mm                       | 450mm                         |
| H0024 | 1500mm      | 250mm                       | 550mm                         |
| H0025 | 1800mm      | 300mm                       | 750mm                         |

### 3. Components & Setup

Included in the package:

- 1 x Pair of Aluminium Curved Non-Folding Loading Ramps
- 2 x Locking Pins
- 1 x Quick Start Guide

#### Setup Instructions:

1. Place both ramps parallel to each other, aligned with the wheelchair or scooter wheels.
2. Ensure the top end lip is resting securely on the step, threshold, or vehicle.
3. Insert safety pins to prevent ramp kick-out.

**NOTE:** Always test with body weight before wheelchair use to confirm placement.

---

### 4. Safe Usage Instructions

- Ensure the surface is dry and flat before using the ramp.
  - Centre the wheelchair/scooter with the ramps.
  - Always have an attendant present when in use.
  - Never exceed the maximum weight capacity or incline height.
- 

### 5. Loading & Unloading

#### 5.1 Loading

1. Align wheelchair wheels with the centre of each ramp track.
2. Push or roll slowly and steadily up the ramp.
3. Avoid sudden movements or twisting wheels.

#### 5.2 Unloading

1. Reverse direction carefully.
2. Guide the chair straight down.
3. Keep an attendant nearby during descent.

**TIP:** Use extra caution during wet weather to reduce risk of slipping.

---

## 6. Maintenance & Care

- Inspect the ramp before each use for cracks or damage.
  - Do not use the ramps if they are damaged.
  - Clean regularly with a cloth and mild detergent.
  - Store indoors in a dry, secure place when not in use.
  - Use proper lifting techniques when handling to avoid injury.
- 

## 7. Safety Warnings

-  DO NOT exceed rated capacity.
  -  DO NOT use with wheels wider than the ramp track.
  -  NEVER stand between ramp tracks.
  -  DO NOT leave the ramp unattended while deployed.
  - ONLY use on firm, flat, stable ground.
  - ALWAYS secure the ramp before use.
  - KEEP bystanders, pets, and children clear during operation.
- 

## 8. Warranty & Support

We stand by the quality and durability of every Heeve® ramp. That's why your ramps are backed by a 12-month manufacturer's warranty for peace of mind.

If you ever have an issue, our friendly support team is here to help.

Product Manufacturer (Heeve®)

Email: [support@heeve.com.au](mailto:support@heeve.com.au)

Web: <https://heeve.com.au/contact/>

Authorised Distributor (Ramp Champ®)

Email: [support@rampchamp.com.au](mailto:support@rampchamp.com.au)

Phone: 1300 913 047

Web: <https://www.rampchamp.com.au/pages/contact-us>

For warranty and support enquiries, contact Ramp Champ®. For engineering or product specification questions, contact Heeve® directly.

---



## Heeve® Aluminium Curved Non-Folding Loading Ramp

1300 913 047

[hello@heeve.com.au](mailto:hello@heeve.com.au)

[www.heeve.com.au](http://www.heeve.com.au)

Authorised Distributor: [www.rampchamp.com.au](http://www.rampchamp.com.au)