Document Version: v1.0

Date of Issue: June 2025

Product Model: Heeve® 700mm Aluminium Cross Over Pedestrian Bridge Ramp

Capacity: 300kg

Brand: Heeve®

Distributor: Ramp Champ®

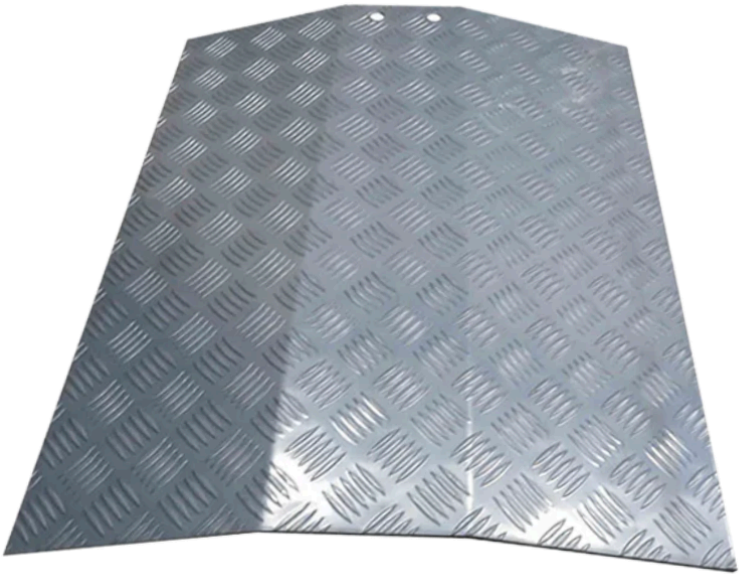


TABLE OF CONTENTS

TABLE OF CONTENTS
1. Introduction
2. Product Specifications
    Gradient Guidelines
3. Components & Setup
    Setup Instructions:
4. Safe Usage Instructions
5. Usage Instructions
6. Maintenance & Care
7. Safety Warnings
8. Warranty & Support

Heeve® 700mm Aluminium Cross Over Pedestrian Bridge Ramp

---

1. Introduction

Thank you for purchasing the Heeve® 700mm Aluminium Cross Over Pedestrian Bridge Ramp. Designed to provide safe and efficient access, this ramp is ideal for crossing thresholds, landings, and small steps with wheelchairs, mobility scooters, and other pedestrian equipment up to 300kg. Please read this manual carefully and follow all instructions to ensure safe, effective use and long-lasting performance.

---

2. Product Specifications

Heeve® 700 Aluminium Cross Over Pedestrian Bridge Ramp	
Ramp Type	Aluminium Cross Over Pedestrian Bridge Ramp
Total Lengths Available	540mm to 980mm
Overall Ramp Width	700mm
Usable Width	700mm
Maximum Capacity	300kg
Ramp Material	Lightweight Aluminium Alloy
Surface	Check place grip
Application	Suitable for manual or powered wheelchairs/scooters and other pedestrian mobility aids

Gradient Guidelines

- Occupied loading height must be maintained at a maximum of 1:6 gradient or 9.5 degrees.
- Unoccupied loading height should be at a maximum of 1:4 gradient or 14.5 degrees.

SKU	Ramp Length	Ramp Width	Channel Height
H0476	540mm	700mm	40mm
H0475	550mm	700mm	50mm
H0477	590mm	700mm	60mm
H0478	780mm	700mm	80mm
H0477	980mm	700mm	100mm

3. Components & Setup

Included in the package:

- 1 x Aluminium Cross Over Pedestrian Bridge Ramp
- 1 x Instruction Manual

Setup Instructions:

1. Place the ramp flat over the threshold or doorway, ensuring both ends sit flush.
2. Press down gently in the middle to confirm it's stable and won't shift.
3. Roll or walk across slowly, maintaining control.
4. Supervise anyone using the ramp to ensure confidence and safety.

NOTE: Always test with body weight before wheelchair use to confirm placement.

---

4. Safe Usage Instructions

- Ensure the surface is dry and flat before using the ramp.
- Centre the wheelchair/scooter with the ramps.
- Never exceed the maximum weight capacity or incline height.

---

5. Usage Instructions

1. Place the ramp securely across the gap, threshold, or step.
2. Ensure both ends of the ramp rest flat and stable on each surface before use.
3. Check that the area is clear of debris, water, or obstacles.
4. Walk, roll, or push mobility equipment (such as wheelchairs or trolleys) slowly and centrally across the ramp.
5. Avoid sudden turns or quick movements while crossing.
6. If possible, have an attendant assist for added safety.

TIP: Use extra caution during wet weather to reduce the risk of slipping.

---

6. Maintenance & Care

- Inspect the ramp before each use for cracks, bends, or damage.
- Clean regularly with a cloth and mild detergent.
- Store indoors in a dry, secure place.

---

7. Safety Warnings

- ⚠️ DO NOT exceed rated capacity.
- ⚠️ DO NOT use with wheels wider than the ramp track.
- ⚠️ NEVER stand between ramp tracks.
- ⚠️ DO NOT leave the ramp unattended while deployed.
- ONLY use on firm, flat, stable ground.
- ALWAYS secure the ramp before use.
- KEEP bystanders, pets, and children clear of the ramp during operation.
- ALWAYS follow manufacturer recommendations for chairs or scooters.

---

8. Warranty & Support

We stand by the quality and durability of every Heeve® ramp. That's why your ramps are backed by a 12-month manufacturer's warranty for peace of mind.

If you ever have an issue, our friendly support team is here to help.

Product Manufacturer (Heeve®)

Email: support@heeve.com.au
Web: https://heeve.com.au/contact/

Authorised Distributor (Ramp Champ®)

Email: support@rampchamp.com.au
Phone: 1300 913 047
Web:https://www.rampchamp.com.au/pages/contact-us

For warranty and support enquiries, contact Ramp Champ®. For engineering or product specification questions, contact Heeve® directly.

Heeve® 700mm Aluminium Cross Over Pedestrian Bridge Ramp

1300 913 047

hello@heeve.com.au

www.heeve.com.au

Authorised Distributor: www.rampchamp.com.au